
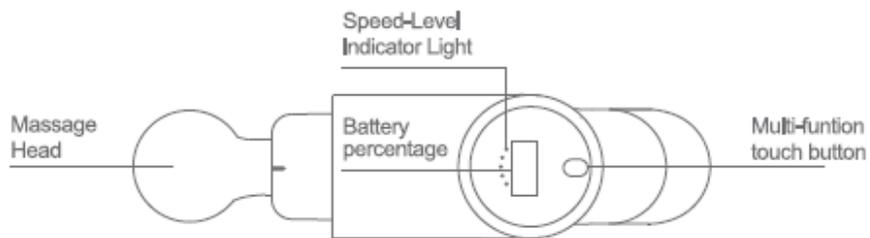
PRODUCT INTRODUCTION

Side View

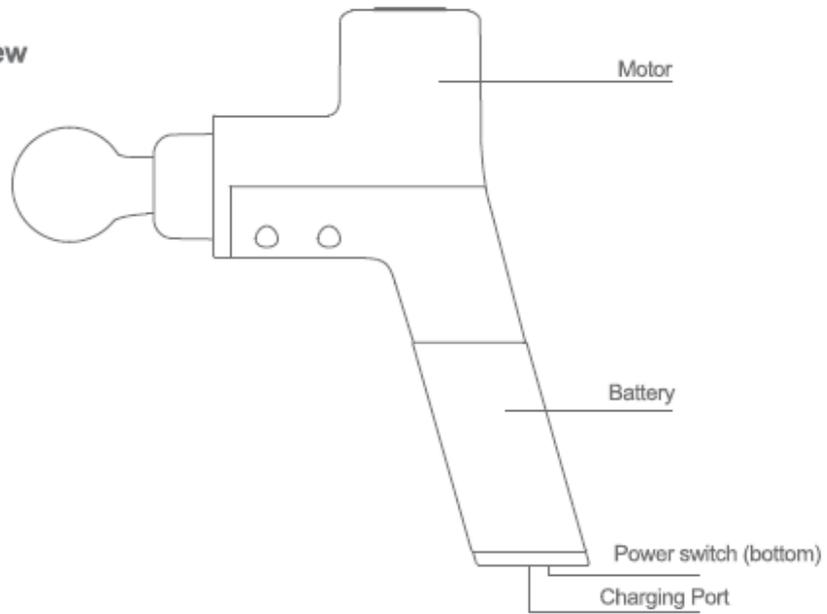Motor

Battery

Power switch(button)

Charging port

Massage head

Speed-level indicator light

Battery percentage

Multi-function touch button

Produkt introduktion

Motor

Batteri

Afbryder (knap)

Opladningsport


Massage hoved

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Batteriprocent

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4 INTERCHANGEABLE HEAD ATTACHMENTS

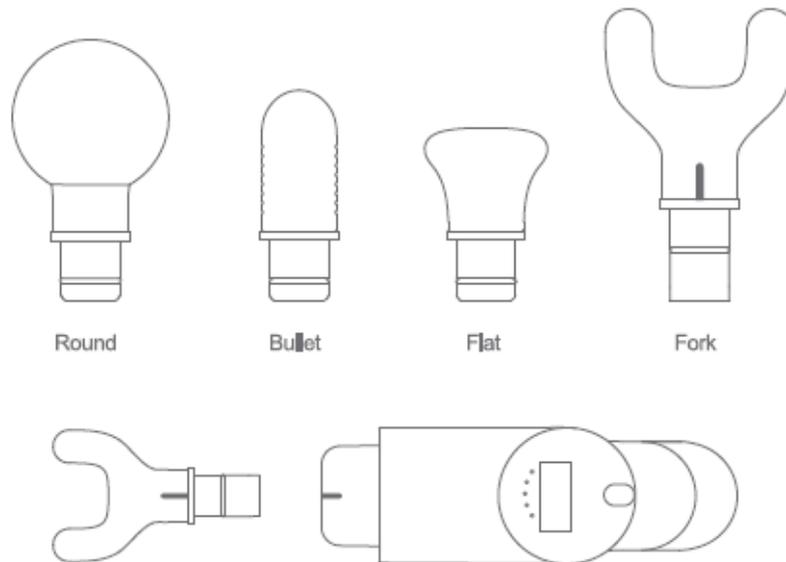4 interchangeable head attachments

Round

Bullet

Flat

Fork

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Rund

Kugle

Flad

Gaffel

Product introduction

Speed-level Indicator Light: level 1 to 5 (left to right)

Battery percentage: 1% to 100%

Power switch (bottom): Main switch to power on/off the machine. Multi-function touch button: start/stop the machine and adjust the speed level.

4 interchangeable head attachments

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4 udskiftelige hovedtilbehør

Assembling: Insert the head into the connector hole. Press and rotate it till the end. Dis-assembling: Pull the head and rotate it simultaneously till it is out.

Specification

Powerful high torque motor with quiet glide technology.

Rechargeable lithium-ion battery, (more than 4 hours playing time)

Ergonomic design and easy to use.

5 speed settings for optimal pressure up to 3,200 RPM.

Cordless and light weight which is portable. (about 1kg)

4 interchangeable head attachments.

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Trådløs og letvægts (ca. 1 kg)

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WARNINGS AND REMINDERS

Keep away from liquids or heat sources. Keep the vents free of dust and debris.

Remove the charger from battery after fully charged.

Please use the device according to the instructions for use.

Check the machine carefully before use.

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PRODUCT PARAMETERS

Product name: Sports Massage Gun	
Model: HC1244	Product weight: 0.86kg
Work voltage: 12V	Product size: 25*16.8cm
Rated power: 120W	Packing list massage gun*1, pouch*1, massage head*4, user
Battery capacity: 2000mAh	Product certification: CE.ROHS.FCC

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Batterikapacitet: 2000 mAh Produktcertificering: CE.ROHS.FCC

IMPORTANT SAFEGUARDS

It is forbidden for minors and the injured to use this machine.

It can be used to recover muscle damage, but please consult with your doctor before use.

It is forbidden to use on wet skin. It can be used across the clothes.

Do not use it above the neck or any hard and boned part of body. Use the product only on the soft tissues of the body.

Do not immerse the machine in water.

Do not use a non-DC 24V charger to charge the machine.

Do not disassemble the machine in any way.

Do not put fingers, hair or other parts of your body get close to the shaft and back of the vibrating head.

Do not touch this part before starting the machine without installing the massage head.

Do not touch this part when start the machine after installing the massage head.

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OPERATING INSTRUCTION

Charging

Fully charge the machine before first use.

Use the USB-C cable for charging.

The power indicator light will flicker quickly when charging.

The LED screen shows the balance battery percentage.

The product cannot be turned on or used when charging.

Please charge the product when less than 5% battery capacity.

The average playing time of the product is about 4 hours.

Operation

Open the power switch at the bottom of the handle. Then the balance battery percentage will be shown on the LED screen at the top of the product.

Short press multi-function button to make the product running. It is on 1st speed level from start and each touch is for one more speed level. There are 5 speed levels in total and the product will stop running when you touch 6 times.

Turn off the product by power switch. The LED screen and indicator light will be off.

Note:

- The product can be turned off by power switch under any speed level.

The product will start intelligent protection and stop running after 10 minutes of continuous use.

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